


FS-3200 Real-time Fluorescent Quantitative PCR Instrument

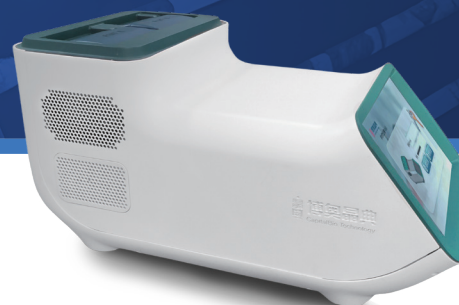

● Product Introduction

FS-3200 Real-time Fluorescent Quantitative PCR Instrument is a medium throughput rapid real-time fluorescent quantitative PCR instrument independently researched and developed by Beijing CapitalBio. Based on the real-time fluorescent quantitative PCR test principle, it is used for qualitative and quantitative test or melting curve test of target nucleic acid (DNA/RNA) in sample. FS-3200 adopts the design of double modules and 16x2 wells, and modules are non-interfering and can be configured independently for operation, maximizing the space utilization. As a medium throughput rapid PCR instrument, FS-3200 can be widely used for clinical test of pathogenic organism, oncogene, reproductive inheritance, drug-related gene and other fields.

● Product Features

01

Multiplex Test

6-color fluorescent channel, it can realize the concurrent testing for multi-index pathogens

02

Simple Operation

The temperature rising rate can be up to 12°C/S

03

Rapid

The test throughput for 1-32 samples can be chosen flexibly


04

Flexible Throughput

The modules can be configured and operated independently, and used efficiently

● Main Performance Index

Optical channel: 6-color fluorescent channel	Temperature accuracy: $\leq 0.1^{\circ}\text{C}$
Sample volume: 1-32	The maximum melting temperature range of melting curve is 35.0°C-95.0°C
Maximal temperature rising rate: $\geq 10^{\circ}\text{C/s}$	The coefficient of variation (CV, %) of fluorescence intensity test repeatability is $\leq 1.5\%$
Maximal temperature dropping rate: $\geq 6.5^{\circ}\text{C/s}$	The coefficient of variation (CV, %) of fluorescence intensity test precision is $\leq 3.5\%$
Module temperature control accuracy: $\leq 0.1^{\circ}\text{C}$	Weight: 12kg
Temperature resolution: $\leq 0.1^{\circ}\text{C}$	Instrument size: 490mm × 265mm × 250mm

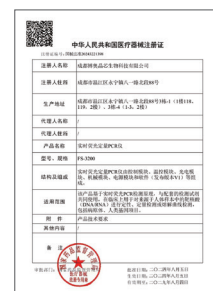

GXZZ 20243221398

CapitalBiotech

Add: 88 Kechuang 6th Street, Building C, Yizhuang Biomedical Park, Beijing Economic-Technological Development Area, Beijing 101111, China
www.capitalbiotechnology.com Tel: +86-10-69002900 E-Mail: globalsales@capitalbiotech.com